



— Pharma Process & Packaging

www.saintyco.com

Vertical Capsule Polisher

C&C 100C Series



CATALOGUE

Vertical Capsule Polisher

C&C100C series Capsule Polisher has achieved double functions of capsule polishing and elevating. Inlet of the machine can be connected with any model of capsule filling machine, and outlet with capsule sorter and metal detector, so that to realize the production line of capsule polishing, elevating, sorting and detecting.



C&C100C

Features:

- Contact parts 316L.
- Up to 600,000 capsules /hr.
- Parts in contact with medicine(main shaft, polishing barrel, feeding port, discharging-port and the channel of sorter) all made of 316L stainless steel
- 360° circle adjustment for inlet and outlet to make production more convenient and save factory room;
- Detachable brush on main shaft can be thoroughly cleaned.
- Automatically sort less dosage, empty, fractured and cap-body separated capsule;
- No blind spot during whole machine cleaning.
- Capsules from the entrance up to the outlet, no retention capsule in polishing barrel.
- Fast-joint system to ensure safe and quick disassembly and assembly;
- cGMP compliant
- CE marked



C&C100CS




C&C100CD-L



C&C100CDS-L

C&C100CS Dimensions:

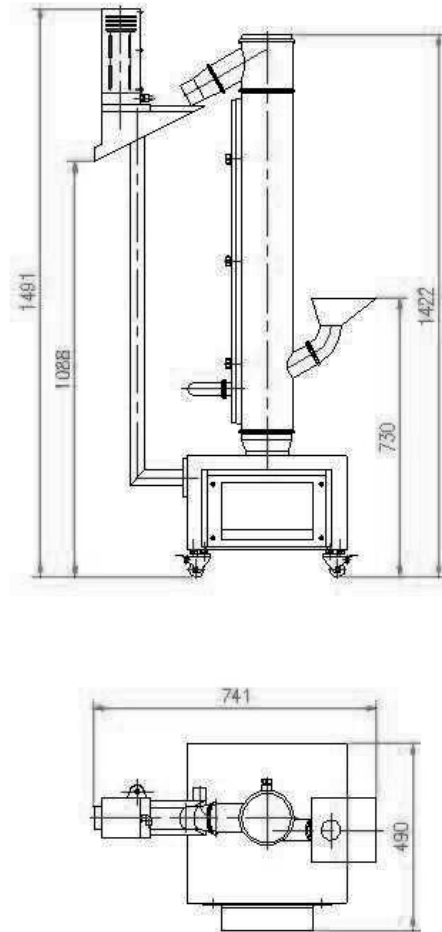

Technical Specifications	C&C100C	C&C100CD-L	C&C100CDS-L	C&C100CDS1-L	C&C100CS	C&C100CS1
Capsule Model No.	00# 0# 1# 2# 3# 4#					
Max Output Capacity (pcs/h)	600,000		300,000	600,000	300,000	600,000
Inlet Height (mm)	730					
Outlet Height (mm)	1260	890	750	800	1050	1100
Power Supply	1PH 220V/110V 50Hz/60Hz 0.2kw					
Compressed Air	0.3m ³ /min 0.3Mpa					
Vacuum	3.0m ³ /min -0.1Mpa					
Dimension (mm)	490 x 510 x 1420	740 x 510 x 14200	1220 x 510 x 1420		740 x 510 x 1500	
Weight (kg)	70	110	120	125	75	80

SaintyCo also offer:

Encapsulation Solutions

Tablet Compression Solutions

Coating Systems

Cartoning Systems

Blister Packing Solutions

